

CMV pp65 polyclonal antibody

CATALOG NUMBER: PAB0003

LOT NUMBER: #

IMMUNOGEN: CMV recombinant antigen pp65 (Rekom Biotech catalog reference RAG0016)

SPECIFICITY: Monospecific, reacts only with CMV*. Non-reactive with normal human serum.

PREPARATION: Isolated from sera of rabbits immunized with highly pure (>98%) recombinant pp65. Anti-CMV pp65 specific antibody was purified by protein G affinity chromatography.

PURITY: >95% (SDS-PAGE)

PRESENTATION: liquid

SOURCE: rabbit

DATE OF BATCH RELEASE: 04/03/2021

BATCH COMPOSITION:

COMPONENTS	COMPOSITION
Storage buffer before lyophilisation	20 mM phosphate buffer pH 7, 0.15 M NaCl and 0.1% sodium azide

LOT SPECIFICATIONS:

- 1. CONCENTRATION (A_{0.1%} = 1.4):** 2.28 mg/ml
- 2. TOTAL QUANTITY PER ALIQUOT:** 1 mg
- 3. TOTAL VOLUME PER ALIQUOT:** 0.438 ml
- 4. STORAGE:** Antibody is shipped with dry ice. Upon arrival, it should be aliquoted to avoid repeated freezing and thawing cycles and stored at -20°C to -80°C. In order to defrost the antibody, maintain the aliquot at 25°C without shaking to avoid aggregation.
- 5. WESTERN BLOT:** suggested titer approx.1:4,500 to 1:5,000

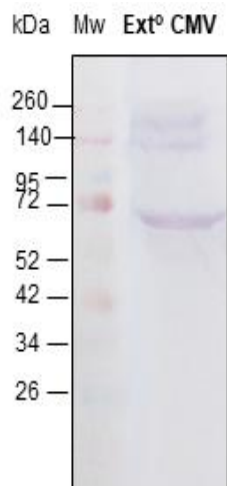


Figure 1. Western blot analysis of the polyclonal Ab to detect the native pp65 in a clarified cell lysate containing a mix of CMV antigens from all stages of the replication life cycle*. Pp65 polyclonal Ab was used in a 1:4,500 dilution over 3 µg of clarified extract.

6. ELISA ANALYSIS:

Antibody titration assay. A titration assay of the antibody was performed with a fixed recombinant immunogen concentration.

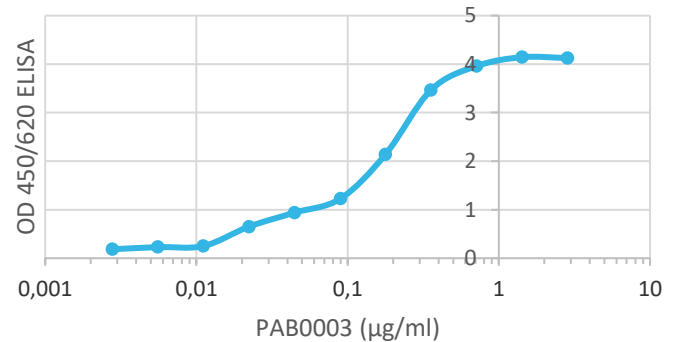


Figure 2. Titration of the polyclonal Ab diluted from 1/800 in a twofold series against constant antigen (1 µg/ml, 19.26 nM). Mean values of OD are used.

The suggested titer for Ab PAB0003 is between 1/51,200 (0.04 µg/ml) and 1/102,400 (0.02 µg/ml).

7. OTHER TECHNIQUES: not tested. Where this product has not been tested for use in a particular technique, this does not necessarily exclude its use in such procedures.

Important Notes: Where this product has not been tested for use in a technique, this does not necessarily exclude its use in such procedures. Suggested working dilutions are given as a guide only. It is recommended that the user titrates.

During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µl or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the containers cap.

FOR RESEARCH AND COMMERCIAL USE *IN VITRO*: not for human *in vivo* or therapeutic use

CMV pp65 polyclonal antibody

Rekom Biotech S.L. – BIC-Granada, Avda. Innovación, 1 – 18016, Granada (Spain) Tel: +34 958 63 70 85
E-mail: info@rekombiotech.com – Web: www.rekombiotech.com – An ISO 9001 Certified Company