

Recombinant antigen ospC for *Borrelia afzelii*

CATALOG NUMBER: RAG0042

LOT NUMBER: #

RECOMBINANT ANTIGEN: the outer surface protein C (ospC) of *Borrelia afzelii* (Lagal *et al.*, 2003).

DESCRIPTION: the *Borrelia afzelii* antigen ospC has been prepared as a recombinant antigen fused to a his-tag in its N-terminus. It is produced from the complete ORF of the gene of the outer surface protein C.

PRESENTATION: liquid protein solution

SOURCE: *Escherichia coli*

MOLECULAR WEIGHT: determined by SDS-PAGE, the protein band is between molecular markers of 25,000 and 35,000 Da, while relative molecular mass calculated from amino acid sequence is 31,320.3 Da.

BATCH COMPOSITION:

COMPONENTS	COMPOSITION
his-ospC	recombinant antigen with a his-tag in its N-terminus
Storage buffer	20 mM phosphate buffer pH 6, 0.3 M NaCl and 0.1% polyoxyethylene (10) tridecylether

QUALITY CONTROL:

1. PROTEIN CONCENTRATION DETERMINED ESPECTROPHOTOMETRICALLY

DO₂₈₀ = 0.132
 A_{0.1%} (=1 g/l) = 0.099
 CONCENTRATION*: 1.34 mg/ml

* The measurement of the protein concentration has been performed with the theoretical extinction coefficient of the recombinant protein obtained from Gill and vonHippel, 1989

2. PURITY CONTROL IN SDS-PAGE: 15%

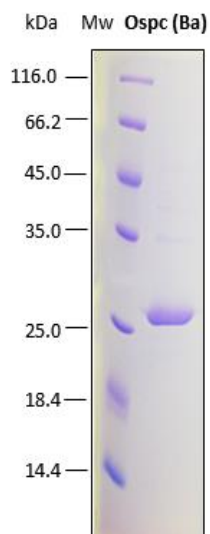


Figure 1. SDS-PAGE analysis (15%) of 5 µl of recombinant ospC. Purity is > 95% as determined by gel electrophoresis. Bands which slightly appear under the main band correspond to this same protein as it is showed in a western blot performed with a his-tag monoclonal antibody.

3. ABSENCE OF PRECIPITATION AFTER A FREEZING AND THAWING CYCLE: ok

LOT SPECIFICATIONS:

- 1. CONCENTRATION:** 1.34 mg/ml
- 2. TOTAL QUANTITY PER ALIQUOT:** 1 mg
- 3. TOTAL VOLUME PER ALIQUOT:** 0.784 ml

4. STORAGE: Protein is shipped with dry ice. Upon arrival, it should be aliquoted in order to avoid repeated freezing and thawing cycles and stored at -20°C to -80°C.

6. APPLICATIONS: not tested. Where this product has not been tested for use in a particular technique this does not necessarily exclude its use in such procedures. Suggested working dilutions are given as a guide only. It is recommended that the user titrates.

7. OBSERVATIONS: Proteins should be maintained frozen at high concentrations. The dilution to be performed for ELISA assays should be made with a small quantity of protein, the same day of the experiment. In order to defrost the protein, maintain the aliquot at 25°C without shaking to avoid aggregation. Prior making test dilutions and after defrost the protein, is recommended to remove possible protein aggregates by centrifuging the stock solution, avoiding alterations in the immobilization of the biomolecule to the solid surface.

RELATED PRODUCTS:

ospC Bb, p41 Ba, p41 Bg, p41 Bb, VlsE Bg, VlsE Ba, ChimLyme Bb.

BIBLIOGRAPHY:

Lagal, V., le Postic, D., Ruzic-Sabljić, E., and G. Baranton. Genetic Diversity among *Borrelia* Strains Determined by Single-Strand Conformation Polymorphism Analysis of the ospC Gene and Its Association with Invasiveness. 2003. *Journal of Clinical Microbiology*, 41:5059–5065

Gill SC, von Hippel PH. Calculation of protein extinction coefficients from amino acid sequence data. *Anal Biochem*. 1989 Nov 1;182(2):319-26.

Important Notes: During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µl or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the containers cap.

Although recombinant antigens are expressed in non-pathogenic *E. coli* and bacterial integrity is destroyed during purification, the antigen preparation should be handled as potentially infectious.

NOT FOR DIAGNOSTIC USE, FOR RESEARCH USE ONLY

Recombinant antigen ospC for *Borrelia afzelii*

Rekom Biotech S.L. – BIC-Granada, Avda. Innovación, 1 – 18016, Granada (Spain) Tel: +34 958 63 70 85
E-mail: info@rekombiotech.com – Web: www.rekombiotech.com – An ISO 9001:2008 Certified Company

Page 2 of 2